


Aspire Lifecomfort Acute Care Mattress



The LIFECOMFORT® ACUTE CARE mattress has three foam layers including a castellated patient surface providing users with immersive, pressure relieving support and air flow for enhanced skin microclimate. The U-shaped base and sidewall structure provides support to assist with patient transfers as well as creates roll resistance at the sleep surface edges. Fitted with superior, hospital-grade polyurethane coated covers designed to withstand the rigors of usage and cleaning requirements in Home Care, Aged Care and up to Acute Care environments.

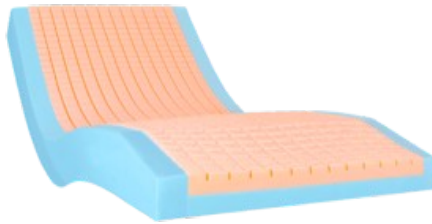

- **Durable Hospital Grade Foam** that is treated for antimicrobial and antibacterial properties
- **Breathable, 4-Way Stretch Waterproof Cover** to reduce shear, friction and maintain skin microclimate
- **Medical Grade Waterfall Cover** is fully-welded over the zipper to maximise infection control
- **Asset Management Ready** with unique batch and serial numbers printed on all covers
- **Carry Handle** - Integrated handles for safe manual handling
- **Hospital Grade** - Suitable for bleach concentrations up to 10,000 ppm
- **4-sided Zipper** - Easily removed for effective cleaning
- **4-way Stretch** - For shear mitigation and friction reduction
- Immersive surface contours and cradles around the user for personalised pressure relief and comfort
- Machine Washable Covers: Up to 95°C temperature.
- Ultra-soft Heel Protection Zone improves pressure relief for vulnerable points. By gently tapering the heel area, pressure is redistributed effectively onto the less vulnerable calf area.
- Castellated Surface collapses in a wave-like movement reducing shear.
- **Foam Layers:**

- IMMERSIVE CASTELLATED SURFACE - Castellated cells offer greater airflow and transitional support
- COMFORT LAYER - Visco-elastic comfort and envelopment
- SAFETY & SUPPORT BASE - Strong U-shaped base for maximum durability and assistive support
- INTEGRATED SIDE WALLS - Strong rigid side layers for transfer efficiency and safety
- **Applications** - Pressure Care: Medium - High
 - Hospital
 - Aged Care
 - Home Care
 - Palliative Care

Warranty Information

- 5 yrs Foam
- 4 yrs Cover

Code	Display Name	Total Length	Total Width	Total Height	Safe Working Load	Unit Weight
NZ-BEM046615	Lifecomfort Acute Care Mattress - Single	1980 mm	880 mm	150 mm	30 -250 kg	14 kg
NZ-BEM046616	Lifecomfort Acute Care Mattress - King Single	2030 mm	1050 mm	150 mm	30 -250 kg	14 kg