


## Aspire Hybrid Comfort Mattress

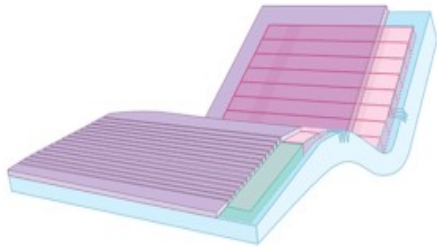

The HYBRID COMFORT is an Acute mattress system offering foam-in-cell technology enabling flexible operation as a high comfort foam static support surface or a powered-up alternating pressure support surface. Specifically re-engineered for increased patient comfort, foam layering now effectively supports patient weight ranging from 30 - 250 kg.

- **Immersive Visco-elastic Memory Foam Surface:** Top layer is now twice as thick to offer significant improvement in user comfort
- **Intelligent Hose Management:** When not in use the air-hose can be stored internally, reducing infection control and OH&S risks
- **Integrated Shear Reduction:** Memory foam layer moves independently of the cell surface reducing shear when bed is adjusted
- **Optimised Heel Protection:** Tapered heel zone redistributes pressure to lower leg, away from the vulnerable heel region. The castellated surface collapses in a wave-like movement reducing shear
- **Innovative pump technology:** Reduces noise and vibration for uninterrupted sleep

**Prescribing Therapists:** Click here to download our [Clinical Assistance Guide](#)

## Warranty Information

- **Warranty Information**
  - 5 Years Foam
  - 4 Years Cover
  - 2 Years Pump/cells

Code	Display Name	Total Length	Total Width	Total Height	Safe Working Load	Unit Weight	Pressure Care Rating
NZ-PTM638005	Aspire Hybrid Comfort Mattress - Single	1980 mm	895 mm	170 mm	30-230 kg	-	Medium Risk
NZ-PTM638006	Aspire Hybrid Comfort Mattress - King Single	1980 mm	1050 mm	170 mm	30-230 kg	-	Medium Risk
NZ-PTM638030	Digital Pump	320 mm	230 mm	130 mm	-	3.5 kg	-